




CONFIDENCE IN LIMITED SPACE

STRAUMANN® NARROW NECK CrossFit®

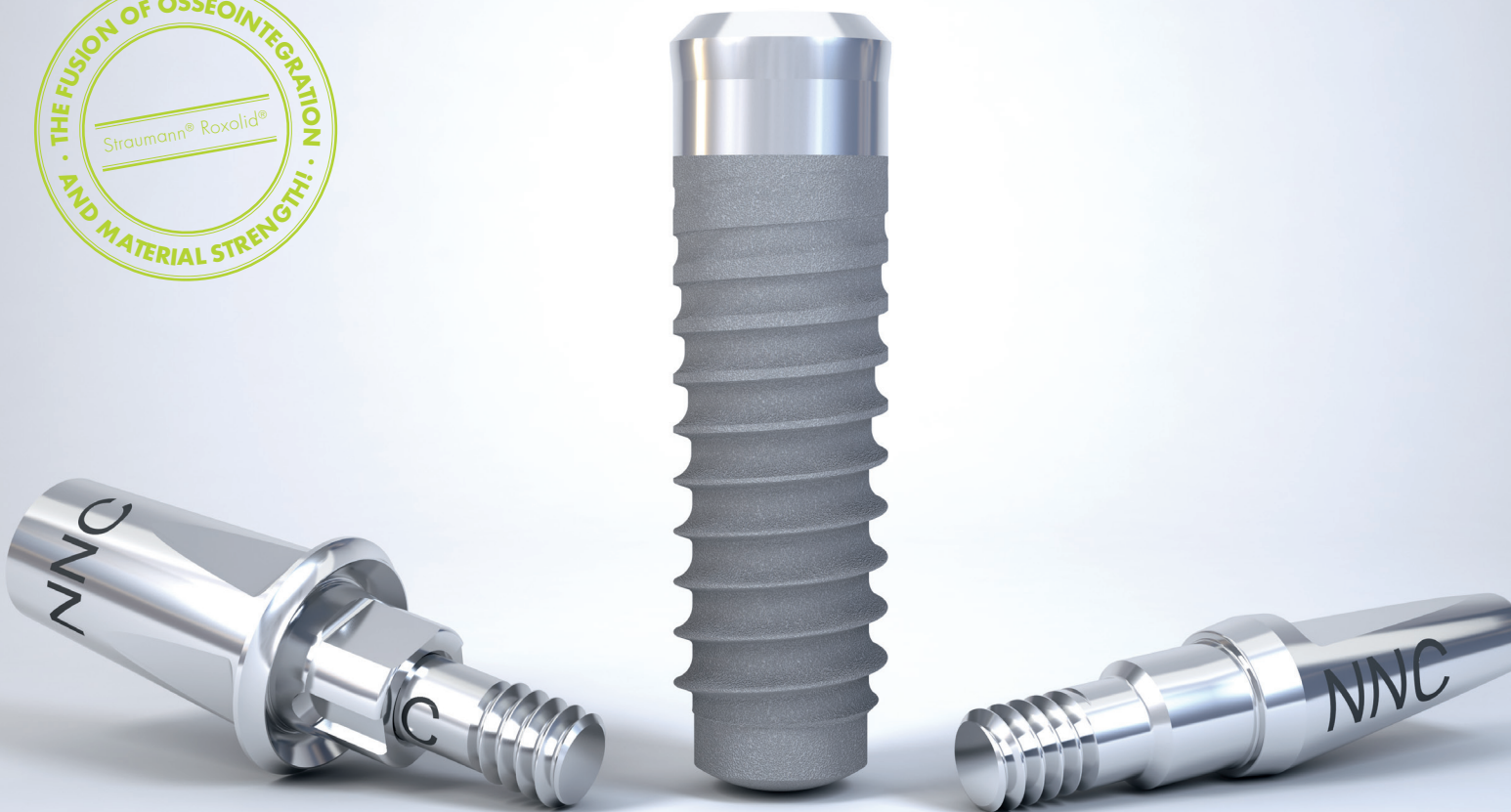

COMMITTED TO
SIMPLY DOING MORE
FOR DENTAL PROFESSIONALSSM


INTRODUCING THE **STRAUMANN® NARROW NECK CROSSFIT® IMPLANT**

THE CONFIDENT SOFT TISSUE LEVEL SOLUTION FOR TREATING LIMITED SPACE.


NNC COMBINES THE LATEST TECHNOLOGIES TO GIVE YOU CONFIDENCE WHEN TREATING NARROW SITES.




Roxolid® material provides high strength and improved confidence when placing narrow diameter implants

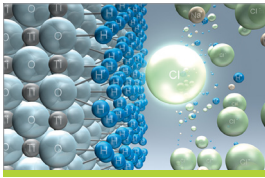

SLActive® implant surface accelerates the osseointegration process¹ and improves predictability in implant treatment

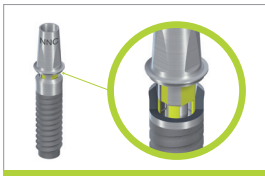

CrossFit® connection offers improved prosthetic flexibility and optimizes abutment insertion²

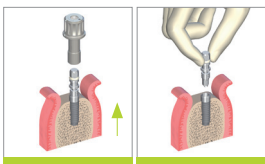

New Transfer Piece simplifies handling during surgical placement

LIMITED SPACE – A CHALLENGE IN YOUR DAILY PRACTICE

Designed for placement in the anterior and premolar region and recommended for use in single-unit, multiple-unit, and fully-edentulous situations³, the NNC can provide excellent treatment outcomes in challenging treatment situations

- Narrow interdental spaces
- Sites with limited bone availability
- Narrow ridges



Clinical example: Lateral incisor in the maxilla

Courtesy of Dr. Roccuzzo



Clinical example: Central or lateral incisors in the lower jaw

Courtesy of Dr. Fragola

¹Compared to SLA® in an animal model

²Compared to Straumann® Narrow Neck

³Small diameter implants not recommended in the molar region



1. CONFIDENCE WHEN PLACING SMALL DIAMETER IMPLANTS


- Roxolid® material
- SLActive® surface

2. WIDE RANGE OF TREATMENT OPTIONS

Comprehensive prosthetic options due to narrow 3.5 mm prosthetic platform with internal connection

3. SIMPLICITY IN DAILY USE


- Integrated soft tissue management through 1.8mm machined collar
- Self-guiding CrossFit® connection at soft tissue level
- New Transfer Piece with snap-fit connection to implant

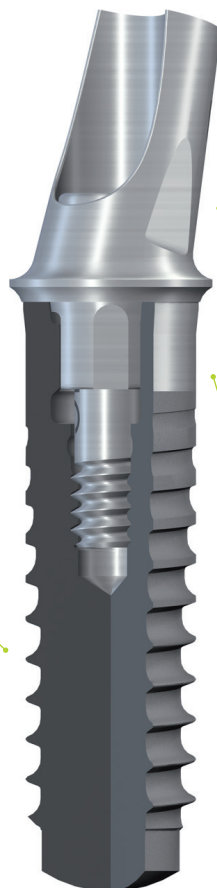

1. CONFIDENCE WHEN PLACING SMALL DIAMETER IMPLANTS

Due to the mechanical strength of the Roxolid® material, it was possible to design an implant with a narrow prosthetic platform and an internal connection while achieving a high Fatigue strength level of 230 N*.

- Straumann® Narrow Neck CrossFit®, 3.3 mm Roxolid® and Straumann® Titanium Abutment
- Straumann® Bone Level NC, 3.3 mm Roxolid® and Straumann® Anatomic Abutment

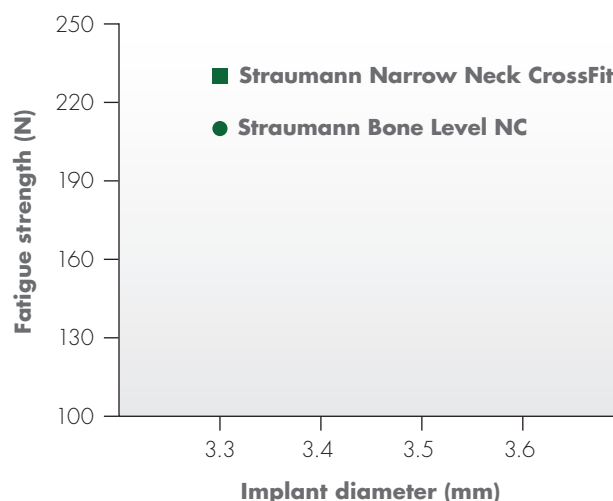
















* Fatigue strength according to ISO 14801 internal tests, data on file (B679A/B567A)

2. WIDE RANGE OF TREATMENT OPTIONS

The Straumann® Narrow Neck CrossFit® Standard Plus Implant System offers a wide range of treatment options:

- Screw-retained and cement-retained restorations
- Choice of implant-level and abutment-level impression workflow

	Screw-retained	Cement-retained			
 Single tooth	 NNC Gold Abutment, crown	 NNC Solid Abutment	 NNC Cementable Abutment, straight	 NNC Cementable Abutment, angled	 NNC Gold Abutment, crown
 Partially edentulous	 NNC Gold Abutment, bridge	 NNC Solid Abutment	 NNC Cementable Abutment, straight	 NNC Cementable Abutment, angled	 NNC Gold Abutment, crown
	Fixed overdentures	Removable overdentures			
 Edentulous	 NNC Gold Abutment, bridge for screw-retained restorations	 LOCATOR®			

3. SIMPLICITY IN DAILY USE

The Straumann® Narrow Neck CrossFit® Standard Plus Implant System is a Soft Tissue Level implant designed to save time and increase efficiency in your practice.

- **Comprehensive prosthetic** portfolio supported by the implant design with integrated soft tissue management.
- Take advantage of the **easy accessibility** to the implant by working at soft tissue level.
- Release the **New Transfer Piece** without using the holding key.



Courtesy of Dr. Roccuzzo

www.straumann.com

International Headquarters

Institut Straumann AG
Peter Merian-Weg 12
CH-4002 Basel, Switzerland
Phone +41 (0)61 965 11 11
Fax +41 (0)61 965 11 01

Straumann USA

Straumann USA, LLC
60 Minuteman Road
Andover, MA 01810
Phone 800/448 8168
978/747 2500
Fax 978/747 2490
www.straumannusa.com

Straumann Canada

Straumann Canada Limited
3115 Harvester Road, 1st Floor
Burlington, ON L7N 3N8
Phone 800/363 4024
905/319 2900
Fax 905/319 2911
www.straumann.ca